

**DAQ box settings**

NI Device:  Input:  Current:  Opto Output:  Voltage:  Sample rate:  Electro Output:

☒ Select default

Measured Impedance:

## E-Stim Parameters

Pulse type:  Optical Stim type:  **29**

Amplitude (mA):  Delay time for e-stim (ms):  **30**

Pulse width (μs):  Light Intensity:  **31**

Interpulse (s):  Optical Pulse width (ms):  **32**

Delay time (s):

Number of pulses:



## Recruitment Curve

Minimum current amplitude (mA):

Maximum current amplitude (mA):

Step size (mA):

Period (seconds):

## Save Options




File name:

☒ Autosave

☒ Minimize Variable Saving

Trial #: 2

Directory: C:\Users\George\Documents\hreflex\_gui\Opto\_Excitatory\_2

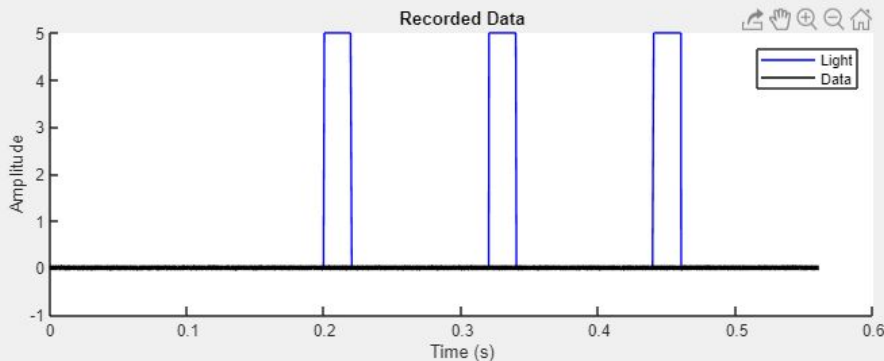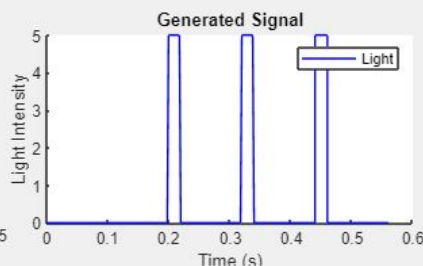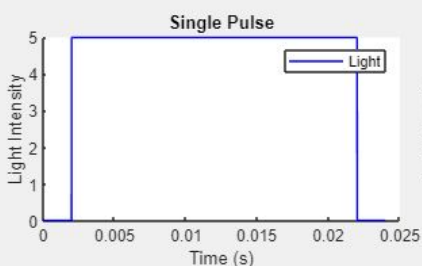

M/H cutoff (ms):

## Explanation of functions (Opto version)

**29** Chooses the optical stimulation type: "Excitatory" outputs only a light pulse, "Inhibitory" outputs a light pulse followed by a brief electrical pulse.

**30** Specifies when the electric pulse is sent after the delivery of the light pulse

**31** Specifies the light intensity: 0 is no light and 5 is maximum light. The LED driver is [turned on at 0] and the option 0 delivers a 5V pulse to the driver, and the option 5 delivers a 0 V to the driver, which turns it on

**32** Specifies the length of the optical pulse in ms